Rotadent[®]

CONTOUR

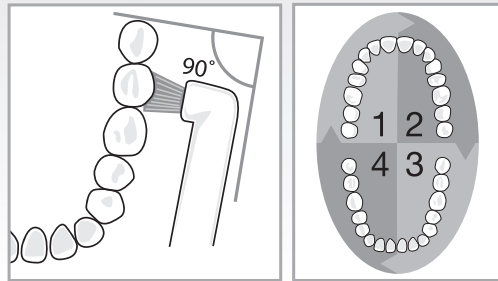
Fig. 1

Fig. 1

Place the brush tip to meet the gum line and tooth surface at a 90 degree angle.

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Fig. 2

Brush each section of your mouth for 30 seconds. A pulse will indicate when it is time to move onto the next section. After the fourth pulse, you will have an additional 30 seconds before the unit turns off.

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Consumers:

To re-order Rotadent brush heads, first contact your Dental Professional, or otherwise visit www.rotadent.com

Includes:

- (1) Powerbrush handle
- (1) One Hollow-Tip brush head
- (1) One Short-Tip brush head
- (1) Charger base with cord
- (1) Travel case

Rotadent[®]

CONTOUR

IMPORTANT

Charge your Rotadent for at least 24 hours prior to use.

This is crucial to the longevity of your Rotadent.

Use a gel-based toothpaste.

This will preserve the useful life of the microfilaments.

Register your product to activate your warranty.

You may register your product at www.rotadent.com

GET STARTED!

1. Charge the unit for 24 hours prior to first use.
2. Use gel-based toothpaste or mouthwash. Paste cleaners will decrease the life and effectiveness of the brush heads.
3. Hold the Rotadent Contour so that the brush tip meets your gum line and tooth surface at a 90 degree angle. (Fig. 1)
4. Cycle through clean and gentle modes by pushing the on/off button. Mode will be indicated by the light.
5. Some individuals may experience a mild tingling sensation, irritation or minor bleeding. If this continues after one week, please consult your dental professional.
6. To clean the brush head and handle, remove the brush head and rinse separately under warm water. Use a damp cloth to wipe down the handle.

For more detailed instructions, clinical references, warranty information and to register your product, visit www.rotadent.com/z_rotadent_how.

DANGER

In order to reduce the risk of electrocution:

- Do not plug the charger into an incorrect voltage source outlet. The Rotadent Contour has been designed to charge with a specific electrical voltage (120 volts @ 60Hz in the United States). Voltage converters and plug adapters do not guarantee voltage compatibility. Never force the plug into an outlet.
- Do not place or store the charger where it can fall into liquid.
- Do not place the charger in water or other liquid.
- Do not reach for a charger that has fallen into water. Unplug immediately.
- Do not use while in the bath or shower.

WARNING

To reduce the risk of burns, electrocution, fire or injury:

- Close supervision is necessary when this product is used by or on children or the disabled.
- Use this product only for its intended use as described in the online manual at www.rotadent.com/z_rotadent_how.
- Do not use attachments other than those provided by the manufacturer.
- Keep cord away from heated surfaces.
- Do not drop the appliance.
- Do not insert any object into any opening of the appliance.
- Do not operate this product if it has been damaged, if it is not working properly, if it has been dropped into water or if the cord or plug is damaged.
- Do not allow your Rotadent Contour to come into contact with any household cleaning products or overspray from any cleaning products used in the area where your Rotadent Contour is stored.

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